

# 3D Perfusion Bioreactor

The transition towards 3D cell culture is revolutionizing traditional cell culture around the world. 3D cell culture provides a more *in vivo*-like environment, thereby allowing the cellular responses from cells cultured in 3D to be more realistic. Furthermore, studies have demonstrated that cells grow significantly better under dynamic culture conditions as a result of the continuous cycling of nutrients, as well as the removal of metabolic wastes. In some instances, the shear force produced by the flowing medium can act as a mechanical stimuli signal that further promotes stem cell differentiation toward certain cell lineages.



As the leader in providing 3D cell culture products and technologies, 3D Biotek is proud to introduce its novel **3D Perfusion Bioreactor**. This unique 3D Bioreactor is a perfect combination of 3D cell culture and dynamic cell culture technologies. The 3D Perfusion Bioreactor consists of multiple independent, autoclavable polycarbonate chambers. The chambers are interchangeable and specially designed to be compatible with 3D Insert™ scaffolds of varying sizes ranging from 96-well to 6-well. Cell culture medias are 100% perfused through the open porous structure of the scaffolds using a pulsatile pump. The entire unit, excluding the pump, is autoclavable and can be used as a single-use bioreactor system.




**Perform Dynamic Cell Culture Using 3D Biotek's Novel 3D Perfusion Bioreactor and Your Choice of 3D Insert™!**

## 3D Insert™ Technology

Porous polymer scaffolds are engineered using 3D Biotek's Proprietary 3D Precision Microfabrication Technology and provide both a 3D cell culture environment, as well as a significantly greater total cell growth area than traditional tissue culture plates. These scaffolds are available in both biodegradable (polycaprolactone, PCL) and non-biodegradable (polystyrene, PS) polymers.




**3D Insert™-PCL Structural Parameters**

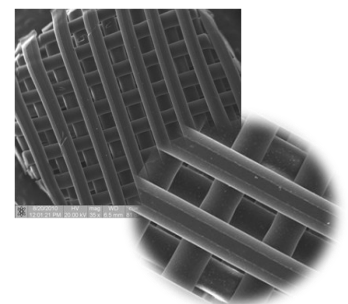
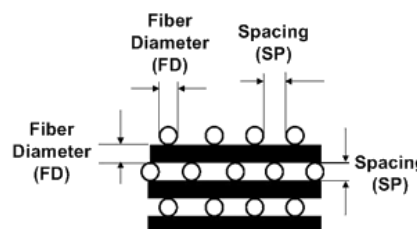
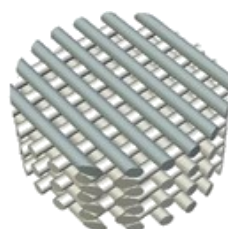


## 3D Perfusion Bioreactor System: *Features & Benefits*

### Chamber/Scaffolds Assembly

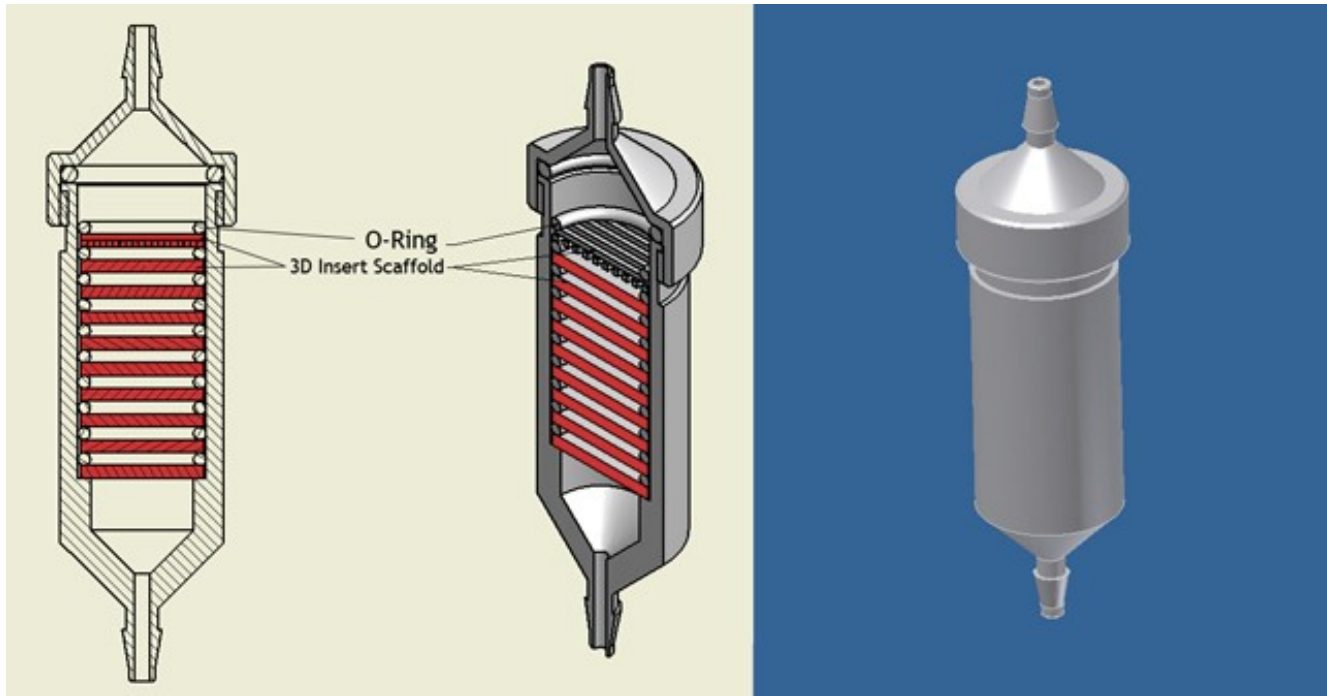

- ◆ The bioreactor consists of four independent, autoclavable polycarbonate chambers that hold up to 10 scaffolds each.
- ◆ The chambers are interchangeable and compatible with 3D Insert™-PCL scaffolds ranging in size from 96-well to 6-well.
- ◆ Within each chamber there is a 1.5 mm distance separating each scaffold. This distance, combined with the offset fiber configuration, ensures the perfusion of media. Slow perfusion of both sterile CO<sub>2</sub> and media is achieved with a low pump speed at the same time.
- ◆ Chambers containing scaffolds can be easily removed and disassembled.
- ◆ The entire unit (except for the pump) is autoclavable and can be used as a single-use bioreactor system.

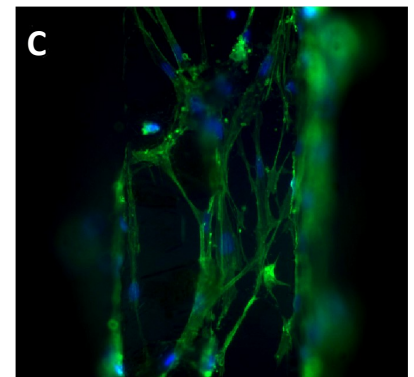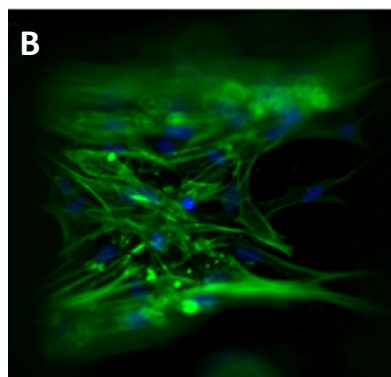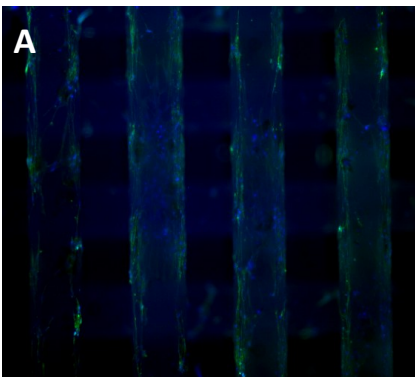

Fluorescent images of hMSC-osteoblasts at Week 1 on 3D PCL 3030 scaffolds (F-actin: green, DAPI: blue, A: 40X, B-C: 200X).